



ECUBE*i7*

Innovation Meets Portability.
Indeed.



Happier Living Everyday

ALPINION MEDICAL SYSTEMS
We are Ultrasound Professionals

E-CUBE i7, the next step in confident diagnosis
with **excellent image performance**.




The innovative portable system

The advantages of a cart-based system made the innovative portable system, E-CUBE i7. To overcome the shortcomings in image quality of a portable system, ALPINION built the E-CUBE i7 around a cart-based system's architecture. Thanks to this powerful and stable architecture, the E-CUBE i7 delivers remarkably excellent and uniform image performance throughout its lifetime.




High-performance with Clinical versatility


■ Elbow joint injection

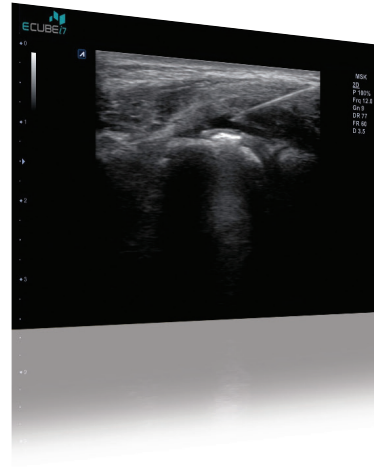

■ Renal perfusion in power Doppler

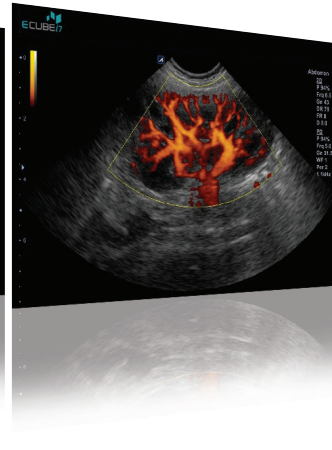

■ Mitral Valve M-mode

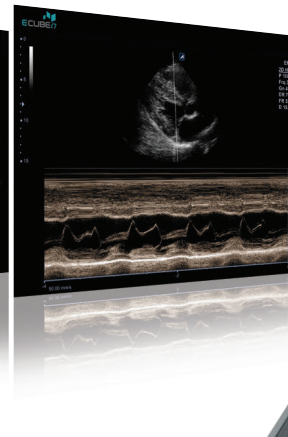

Imaging Technologies


- Advanced Spatial Compounding Imaging (SCI)
- Harmonic Imaging (FTHI/PITHI)
- Speckle Reduction Imaging (SRI/FullSRI™)
- Tissue Doppler Imaging

Software Technologies


- Needle Vision™
- Xpeed™ (Auto Image Optimization)
- Auto IMT measurement
- CW
- Panoramic Imaging



■ Tissue Doppler Imaging

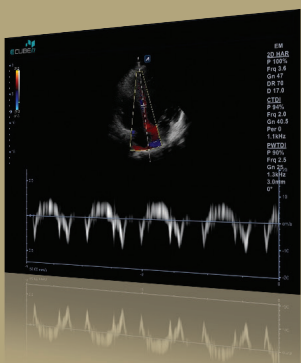


■ Thyroid nodule in Live Dual mode

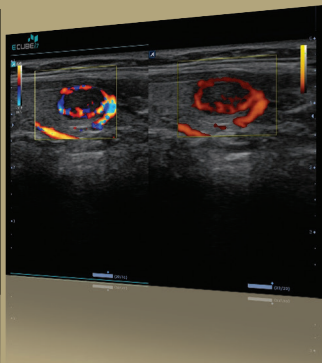


■ Fetal echo with FullSRI and FTHI



■ Umbilical Artery Doppler

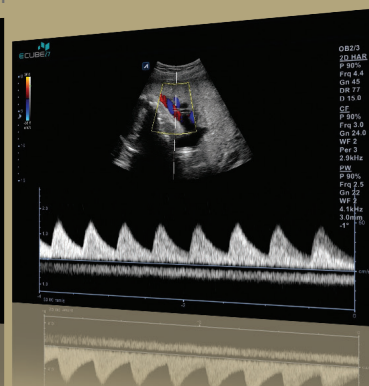

■ Uterus with FullSRI and SCI

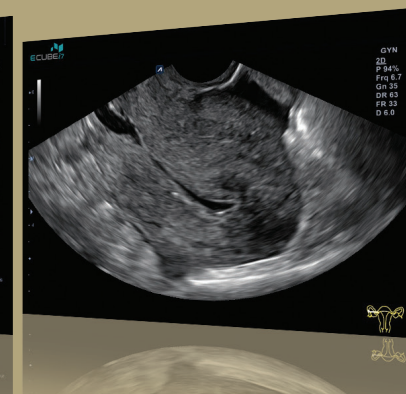


Transducer Guide

Convex	
② C5-8NT	⑥ C1-6T
Linear	
① IO8-17T	⑤ L3-12T
Phased	
③ SP1-5T	④ SP3-8T
Endo-vaginal type	
⑦ EC3-10T	⑧ EV3-10T

The E-CUBE i7 fulfills your long awaited need for **high quality imaging**

Crystal Signature™ : High class single crystal transducer technology


Single crystal is a new kind of piezoelectric material which shows superior image quality. Crystal Signature™, ALPINION's world leading single crystal technology, provides higher sensitivity, better penetration, and less heat degradation than traditional PZT materials. With its single crystal phased array transducers (SP1-5T, SP3-8T), the E-CUBE i7 fully supports various diagnosis including adult echo, pediatric echo, and pediatric abdomen.

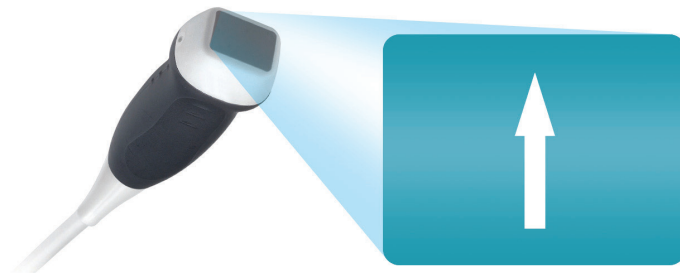


Single crystal material exhibits fewer poling variation




- Increased Bandwidth
- Better Signal to Noise Ratio
- Improved axial resolution and penetration

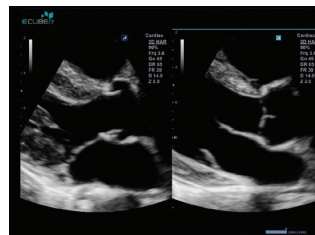


Enhanced contrast resolution

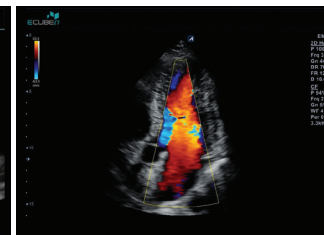

Excellent spatial sensitivity

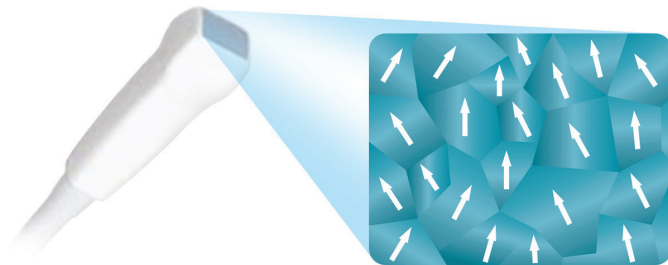


The traditional PZT technology - The variations of polarization in PZT affect its piezoelectric properties and signal to noise ratio(SNR)

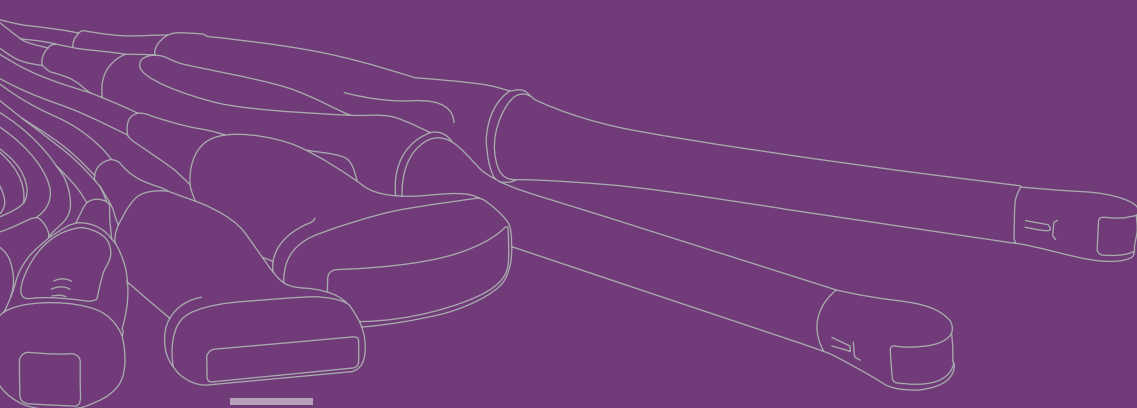

High-frequency enhanced linear transducer


The E-CUBE i7 features a high-frequency linear transducer, the hockey stick-shaped IO8-17T. The IO8-17T is ideal for imaging for superficial, nerve blocks, and vascular procedures. IO8-17T's small footprint allows easier access to tight spots and its excellent resolution in the near field provides extremely detailed anatomical information for your patients.

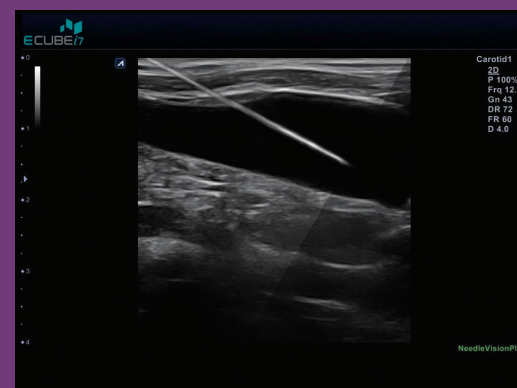

Intra Jugular Vein injection with Needle Vision™ Plus


Portable, **time to slide.**

Experience the slide.
Enjoy the comfort.
Believe in the innovation.

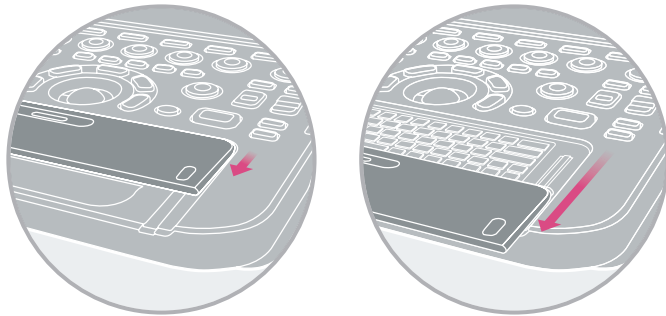

The E-CUBE i7's sliding keyboard cover provides a solution for sonographers who suffer the inevitable wrist pain from the long term operation of imaging platforms. By making small changes, the E-CUBE i7 grants huge benefits. The location of keyboard is closer to user, allowing better support for the arms and wrists, thus reducing wrist pain from unnatural and awkward postures when typing.

The innovative sliding cover supports the wrist at the right angle for a healthy wrist, improves the sonographer's comfort when operating the trackball, and prevents the keyboard from being polluted by dust, gel, fluid, and blood.

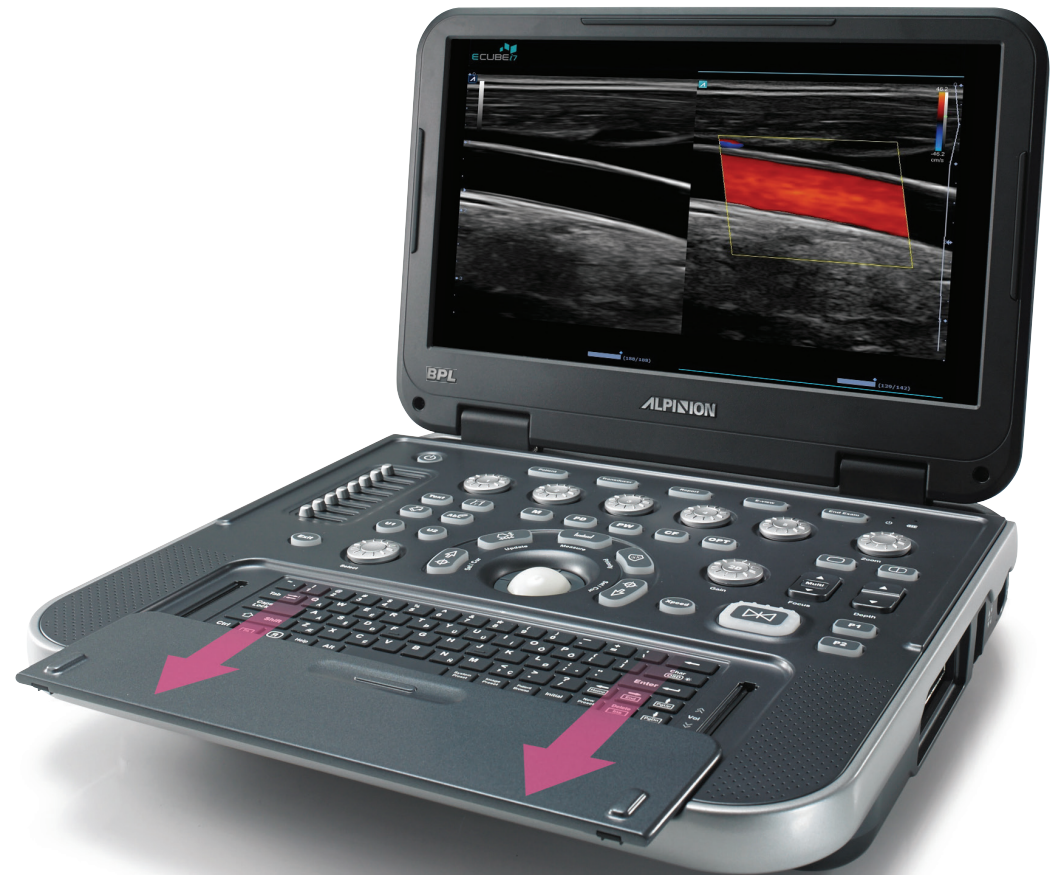


The innovative E-CUBE i7 rewrites the history of ergonomic design for portable ultrasound systems.

High-frequency enhanced linear transducer

Ergonomics you have never seen before.
The E-CUBE i7, the world first portable ultrasound system with a sliding keyboard cover

The E-CUBE i7 features a high-frequency linear transducer, the hockey stick-shaped IO8-17T. The IO8-17T is ideal for imaging for superficial, nerve blocks, and vascular procedures. IO8-17T's small footprint allows easier access to tight spots and its excellent resolution in the near field provides extremely detailed anatomical information for your patients.

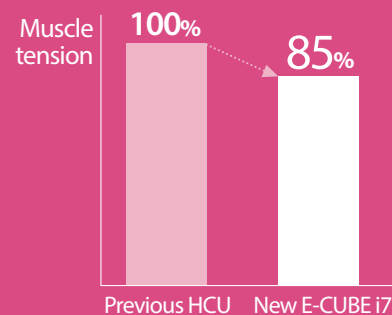

「This study was conducted
by KOREA University Guro Hospital」

Specialized work postures involving the repetitive movement of the joints can lead to an increase in the frequency and severity of repetitive musculoskeletal injuries. Sonographers are at the risk as a result of the posture of their hands, wrists, elbows, shoulders, neck, and upper back.



Previous HCU



E-CUBE i7



Previous HCU



E-CUBE i7


Reinvents **point-of-care**.
Meets the need for **on-the-spot**.
Extends diagnostic **dynamics**.

Workflow




The E-CUBE i7's intuitive workflow helps the user conduct every diagnosis easily and efficiently, improving patient throughput. ALPINION's Xpeed[™] one-touch optimization enhances scanning efficiency by automatically optimizing images in B-mode, Color Flow, and Doppler modes with just the touch of a single button.

Speed

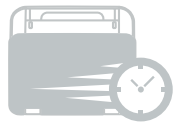

The embedded SSD not only enhances the E-CUBE i7's operation speed, but also its stability. From boot-up to initializing a scan, the time is less than 60 seconds.

Durability

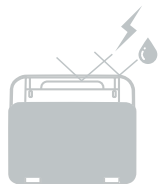

Constructed from solid aluminum, the E-CUBE i7 is built for durability in any situation. The E-CUBE i7's control panel sits on a wide rubber plate that protects its electric boards and inner components from pollutant such as dust, gel, fluid, and blood.



Innovation Meets Portability. Indeed.





ALPINION MEDICAL SYSTEMS

Verdi Tower, 72, Digital-ro 26-gil, Guro-gu, Seoul, 152-848, Korea

Homepage www.alpinion.com

E-mail international@alpinion.com

TEL +82-2-3282-0900

FAX +82-2-851-5591

Copyright©2015 ALPINION MEDICAL SYSTEMS CO. LTD.
All rights reserved.

BPL MEDICAL TECHNOLOGIES PRIVATE LIMITED

11th KM, Bannerghatta Road, Arakere, Bangalore - 560076

Phone +91 80 26484388/ 2648 4350

Toll Free 1800-425-2355

Website www.bplmedicaltechnologies.com

For Enquiries sales.medical@bpl.in

CIN: U33110KA2012PTC067282

Copyright©2015 BPL MEDICAL TECHNOLOGIES PVT.LTD
All rights reserved.

The contents of this catalogue may be partly changed without the prior notice to customer, due to performance enhancement
BPL Medical Technologies is committed to abide by the PC-PNDT law and shall ensure that all regulations applicable to portable Ultrasound devices under the prevailing PC-PNDT law shall be strictly adhered to in every state of the country.
Available for sale in India & Nepal Only.

ECUBEi7:FC:10:15